ERGONOVA

I FEEL GOOD!

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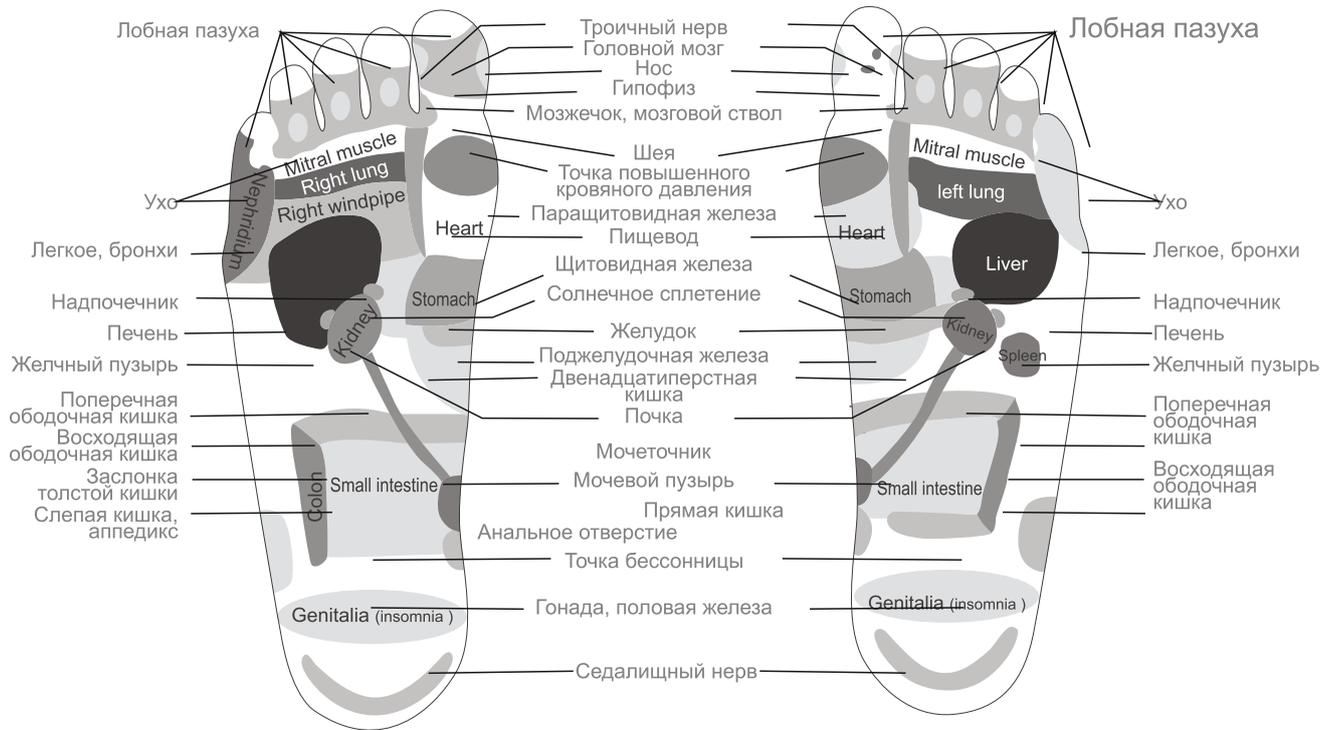
TECHNICAL DATA

Description	FOOT MASSAGER
Model	MaxiGO
Rated Voltage	<input type="checkbox"/> AC 110–120V 60Hz <input type="checkbox"/> AC 220–240V 50Hz
Rated Power	50W
Rated Operating Time	15min
Security Level	<input type="checkbox"/>



ATTENTION

1. This machine is unsuitable for heart disease, patients with blood, broken foot inflammation. Moss feet, all pregnant women and all patients.
2. Please use the product in the rated voltage.
3. Please control the massage time well. That would be better if 30min for each time.
4. It can promote the blood circulation and detox effect, if you drink a glass of water before using the product.
5. Do not use the power cord to drag the foot massage, so as not to cause the power cord injury
6. Don't put this machine in bathroom or moist place, and don't use on the cover with carpeted.
7. Please do not put the product into the water or spray it directly when you clean the product.
8. Please do not reform, disassembly and repair the product if you are not professional worker.
9. Please do not stand in the product directly or put the heavy goods into the machine.
10. The product belongs to family appliance which can not instead of the medical appliance. Please stop using it ,if you feel uncomfortable.
11. The electrical line should be changed by the manufacture, repair department or the similar professional workers ,if the electrical line was broken.
12. Please do not use the product if accessories damaged or exposed internal parts..
13. Children, the disabled or mentally handicapped adults must have the assistance be used.
14. Children and the disabled can not played or tread on the machine.
15. Please use the product carefully if you are insensitive for the heating.
16. Please use the product in the flat ground.
17. Please stop using when the relative component was broken.

SAFETY USAGE


1. Please do not pull the plug when your hand is wet.
2. Please do not make several appliances in the same socket, in order to avoid fire.
3. Please do not use the broken electrical line or socket.
4. Please cut off the power supply to avoid the electrical line and socket from breaking.
5. Please pull up the plug if you will not use it in the long time.
6. Please do not drag or tread the electrical line.
7. Please make sure the plug connects well with the socket. And please do not use it when the electrical line, electrical plug and socket are broken.
8. Please do not use the product when it occurs thunder and lightning, also pull up the plug.
9. The product can not be put in an unstable place. Because it will fall easily due to shaking or external causes.
10. Please do not use essential oils, painting and insecticides etc.
11. Don't use it if your leg or foot is injured, follow your doctor.
12. Store in the reach of children, Prevent children from playing with it.
13. If an exception occurs when using (like: Burning smell), Must stop using and unplug the plug.
14. Do not use near fire and flammable.
15. Do not use for a long time in the sun exposure.
16. Before using, do not wear jewelry in the feet, Do not go to sleep while using.
17. Unwell or drunk, do not use.
18. Do not use it on top of electric blankets or heaters, so as not to cause a fire.

PRODUCT STRUCTURE

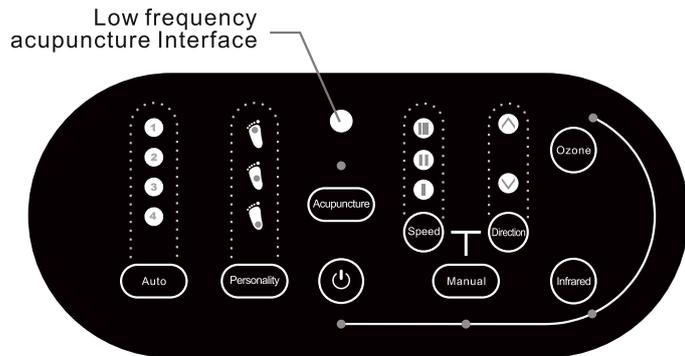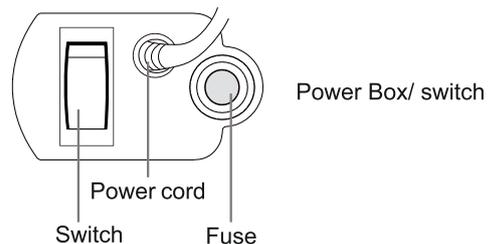
Operating Panel

OPERATION

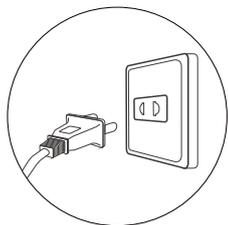
1

Plug in the power.

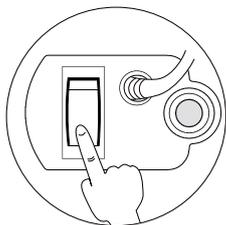
2

Start the switch.

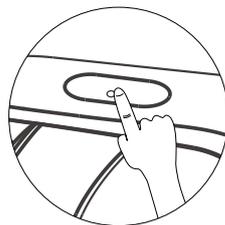3

Put the feet on the location of massage. Press the key on the operation panel to open, select the desired function.

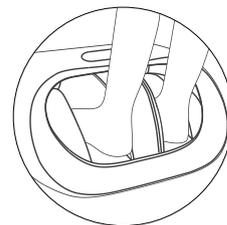
4

Enjoy the massage fun.

The use of low frequency acupuncture

- Put the cable and the butterfly paster buckle together.
- Put the connector into the corresponding jack panel
- Strip off the film on the paster, affixed to the proper location of body.(please clean the skin firstly)
- Start the acupuncture key, which have 3 levels, weak (green light),medium (yellow light),strong (red light).Please start acupuncture under the proper level.

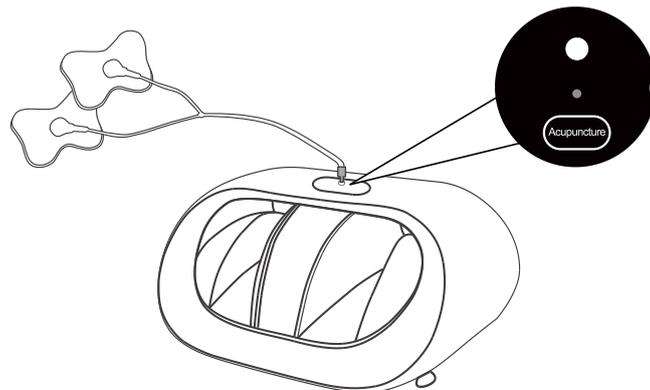The interface cable into the jack corresponding to the operation panel.

FUNCTION

Kneading Massage

Two groups with four clamping piece of kneading massage to stimulate the foot acupoint, promote the blood circulation, improve the sleeping quality, eliminate the fatigue and improve health.

Roller Massager

The massage of the sole roller and plush copper to promote the blood circulation and make the free venation.

Infrared heating

The under of the machine equips with the device of the infrared generator which can produce many kinds of effect. For example, it can promote the blood circulation, improve the metabolism and get the health care.

Low frequency acupuncture


Low frequency acupuncture use the pulse current to setup many kinds of pre-setting program which massage the body effectively. The machine can also use the special position, and remit the fatigue, ache and cervical disease etc. Meanwhile, it is easy to use, and it has health care function for the muscle.

Ozone Sterilization


Built-in the Ozone generator to produce the Ozone to sterilize. Meanwhile, it can clean the cloth cover and foot effectively. It is healthy, environmental and safety.

INSTRUCTION

1.  Power key connect the power and turn on the switch on the power box. When power switch on the operator panel button lights flashing, the machine into standby mode; press the button, power switch button indicator is light, other features on the control panel effectively .

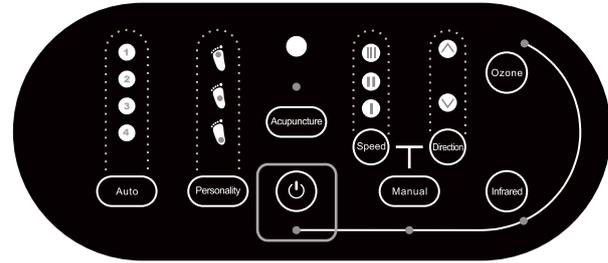
2.  Automatic program control key. There are 4 automatic programs. Each press transfer an automatic program, and corresponding LED will light. When the machine is working in other modes, click this key to enter the automatic mode.

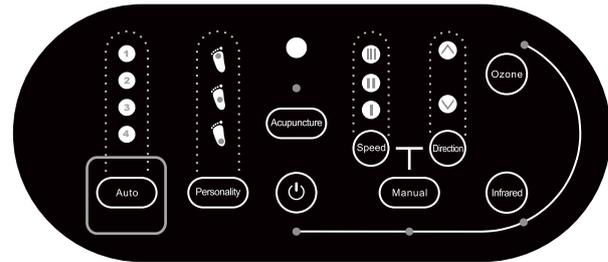
INSTRUCTION

3、**Personality** The key is personality program which has 3 programs. When the machine starts, the key can exit "automatic program" and "manual setting", and access the massage of personality program. It will be another program, when you press the key each time. The related indicator light will be turned on.

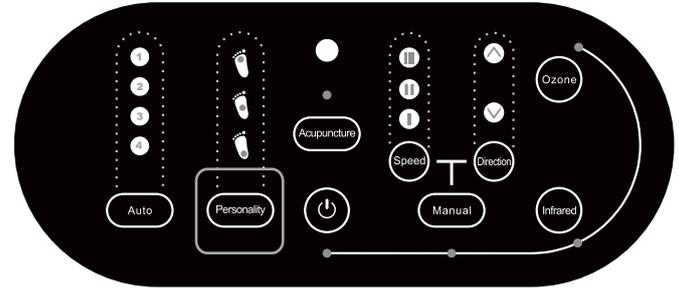
4、**Manual** The manual setting key can control "speed" and "direction". The product has 3 speeds which are low, middle and fast. The indicator light of each speed key has different colors. The direction key has two levels that are clockwise and counter clockwise revolution.

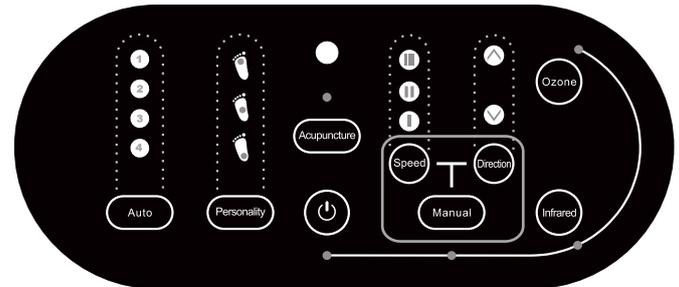
INSTRUCTION

5、 The key has Ozone function. After the machine works, the key open the Ozone function, and the indicator light will be turned on. Meanwhile, the Ozone generator starts working. The function can release the Ozone and sterilize. It will be closed when you press again.

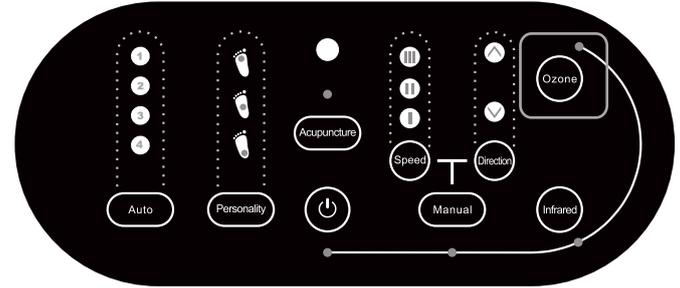
6、 The key has infrared tepidity function. After the machine works, the key can open the infrared tepidity function, and the related indicator light will be turned on. Meanwhile, the internal infrared device will be started to work. The infrared tepidity function will close when you press the key again.

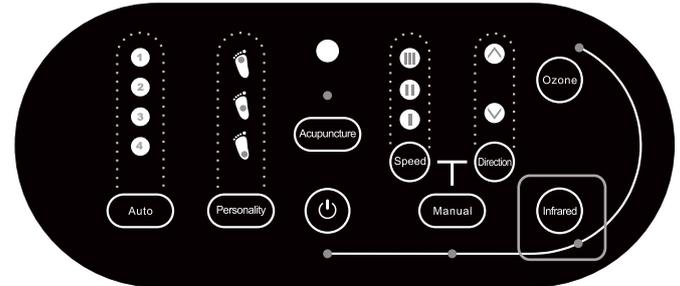
MAINTENANCE

1. Before clean the machine, you should put out the plug.
2. Please don't use the benzene, diluents or other corrosive to clean the machine.
3. If you don't use the machine for a long time, please put it in the clean and dry environment.
4. please do not make it near the heater and avoid exposing to sunlight
5. Make sure to keep in safe when you use the machine, please don't break, press, bump or other heavy strike.

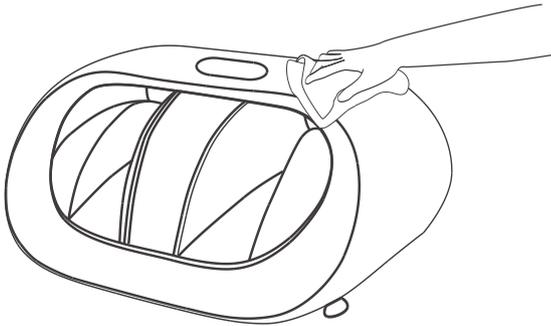○ Please use the soft towel to clean the PU cover

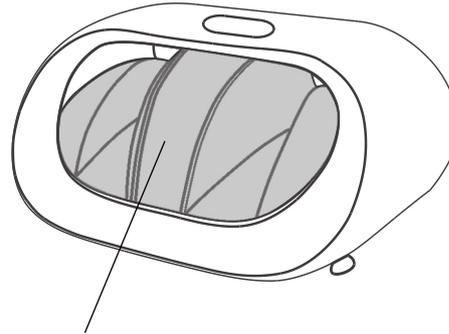○ This is removable and washable fabric cover. when you clean, please open the zipper and take it down, when the fabric cover is dry, you can put it back.


PACKING LIST

1. Foot massager 1PC
2. User manual 1PC

TROUBLESHOOTING

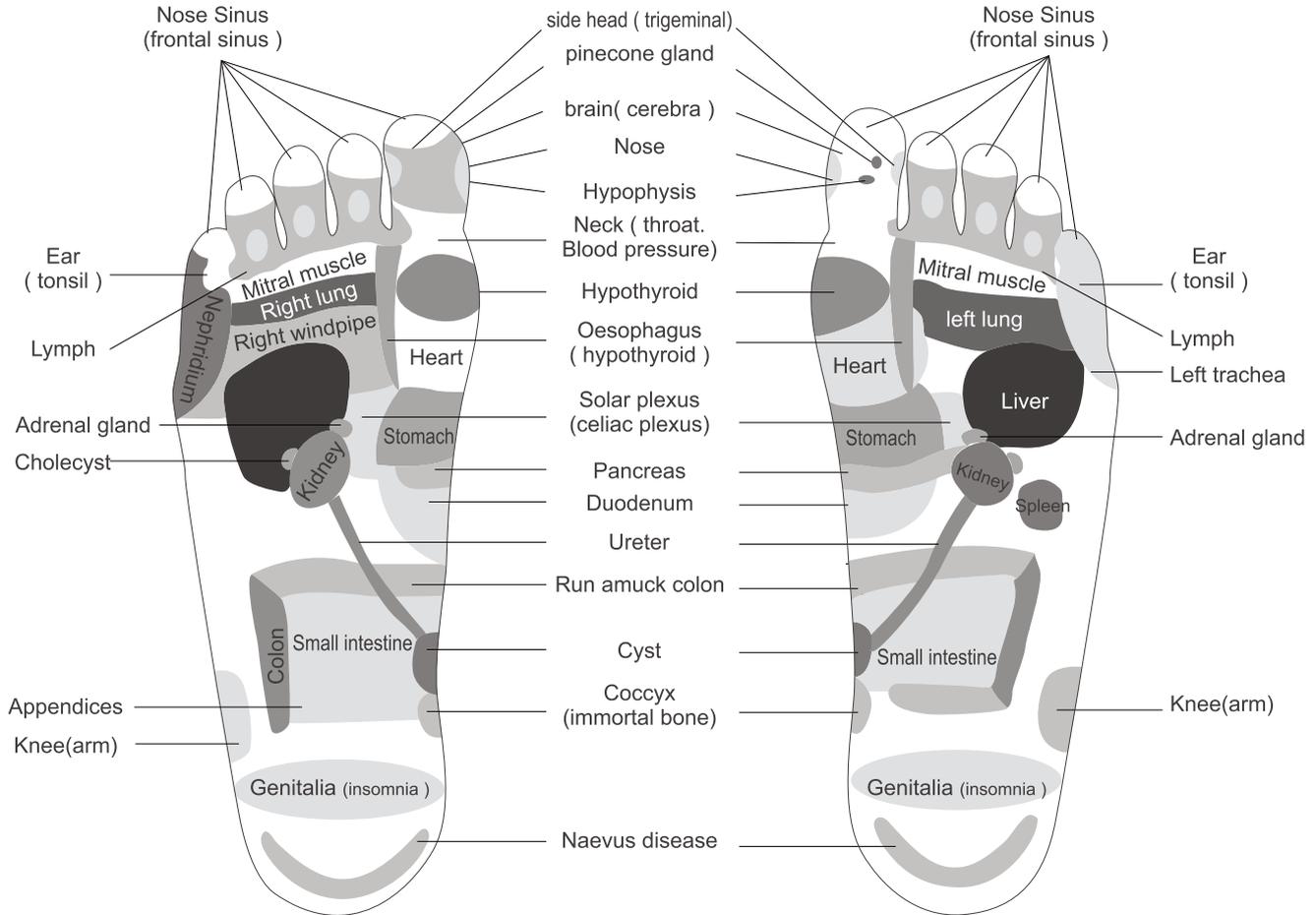TROUBLES	POSSIBLE REASON	SOLUTION
No operation	Switch is in off mode.	Turn on the switch.
	Plug is not correctly inserted into the socket.	Connected to the power.
Stop suddenly	The pressure is too strong.	It will work after reducing pressure.
	In operation for too long. It stops automatically due to the thermal protection function.	Disconnect the power and let it cool for at least 30mins.

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Авторское право © 2012г.